


# SEER™ SALINE ENHANCED ELECTROSURGICAL RESECTION

The SEER resector gives surgeons the flexibility to dissect solid organ parenchyma quickly with reliable hemostasis. Utilizing combined saline and RF current from a standard electrosurgical generator, the SEER provides a zone of coagulation lateral to the incision to minimize blood loss and allow for faster dissection to reduce vascular occlusion times. The wide surface area of the device allows the surgeon to tamponade larger vessels while simultaneously delivering saline and RF current. SEER is available in lengths suitable for open or laparoscopic procedures and fits through 1 cm ports.

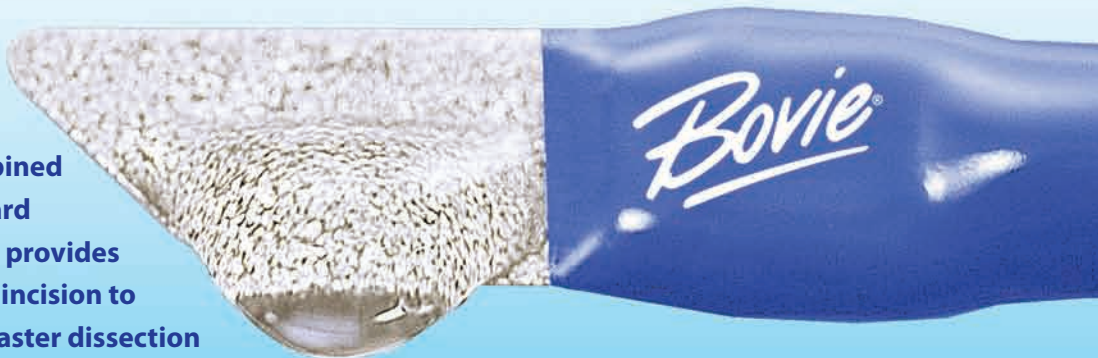

## WHY USE SEER?

- Made with durable materials
- More robust insulation
- Visual and tactile feedback
- Ability to sculpt around lesion
- Economical/Reimbursement Advantage



| PRODUCT NAME      | PRODUCT CODE |
|-------------------|--------------|
| SEER              | SR145L       |
| SEER Laparoscopic | SR326L       |

- SEER delivers saline much closer to the tip
- SEER delivers saline throughout the tip for better saline distribution
- Can be used on any parenchymal tissue
- Large tip allows for tamponading larger vessels
- Optimal temperature distribution
- Multifunctional tip geometry
- Flow control on handpiece
- Potentially shorter vascular occlusion times



- Saline Assisted Cut and Coagulation:
  - Up to 1 cm wide Coagulation Zone
  - 1/2 cm deep cut
- Minimal Charring

